

R REFRIGERATORS



HYC-260L/300L

Model: HYC-260L/300L , Upright Type
Effective Capacity: 260/300Liter
Temperature range: +2°C ~ +10°C



Features:


- Upright type.
- Temperature ranging from +2 to +10°C, can be set freely, controlling precision 1C, display accuracy 0.1C.
- Combined with two exact sensors and one sensors to defrost.
- 2-layer glass door, and inert gas inside.
- Room temperature can be ranging from 0°C to 32°C.
- Forced air circulation system with intelligent control.
- Audible and visual alarm: High or low temperature alarm, Door ajar alarm, system failure alarm.
- 4units Caster are mounted under the bottom.
- Interior fluorescent lighting.
- 5units shelves made of quality steel wire.
- Exterior made from high quality structural plate, interior made from ABS material.
- China good quality compressor (1 unit) and fan motor.
- Optional : temperature printer.
- Refrigerant as R134a, CFC free.
- Power supply: 220V/50Hz 1 phase, can change as 220V/60Hz or 110V/50Hz or 110V/60Hz.
- Certificate: ISO9001, ISO14001,ISO13485.

Model	Capacity	External size (mm)	Internal size (mm)	Input power	Weight (N./G.T)
HYC-260L	260 Liters	D590xW603xH1556	D445xW525xH1070	215	77Kg/105Kg
HYC-300L	300 Liters	D590xW603xH1728	D445xW530xH1265	215	86Kg/110Kg



RXC-240L/280L

Model: RXC-240L/280L,Upright Type
Effective Capacity: 240/280Liters
Temperature range: 4±1°C.

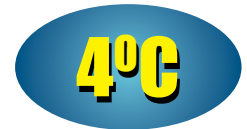

Features:

- Temp. range should have 4±1°C, with temp. printer standard.
- Vertical cabinet (upright).
- Polyurethane insulation should be minimum 80-95 mm.
- Microprocessor-based temperature control system.
- Forced air circulation system.
- Automatic temperature control.
- Large screen LCD display of temperature, display accuracy of 0.1°C.
- Digital display of upper and lower temperature, two sensors. Perfect Audible/visual alarm system: High/low temperature alarm, system failure alarm.
- Exterior made from high quality steel plate.
- 2-ply transparent door with heat preservation.
- Inner doors with 5 units wire shelves.
- Interior fluorescent light.
- China good quality compressor (1 unit) and fan motor.
- The high-effective air-cooled condenser and the fin-type evaporator.
- Refrigerant as R134a, CFC free.
- Power supply: 220V/50Hz 1 phase, can change as 220V/60Hz or 110V/50Hz or 110V/60Hz.
- Certificate: ISO9001, ISO14001, ISO13485.

Model	Capacity	External size (mm)	Internal size (mm)	Input power	Weight (N./G.T)
RXC-240L	240 Liters	D615xW605xH1615	D445xW525xH1070	500	96Kg/105Kg
RXC-280L	280 Liters	D615xW605xH1780	D445xW530xH1265	500	105Kg/110Kg



HYC-950L

Models: HYC-950L, Upright Type
Effective Capacity: 950 Liters
Temperature: +2 to +10°C



Features:

Upright type • Temperature ranging from +2 to +10°C, controlling precision 1°C, display accuracy 0.1°C, with temperature printer (Optional in 260L/300L) • Outer is steel board, inside is stainless steel • Fixed with two exact sensors and one sensors to defrost • Door with heater to prevent ice leak • Audible and visual alarm: High or low temperature alarm, Power failure alarm, Low battery • Two layer glass door, inert gas inside, alarm, Door ajar alarm, Filter blocking alarm, system failure alarm • 4 units Caster are mounted under the bottom • Interior fluorescent lighting • 8 units shelves made of quality steel wire soaked with plastic • Refrigerant as R134a, CFC free • Imported famous brand Danfoss compressor (lunit) and EBM fan • Power supply: 220V/50Hz 1 phase, can change as 220V/60Hz or 110V/50Hz or 110V/60Hz • Certificate: ISO9001, ISO14001, ISO1348S.

Model	Capacity	External size (mm)	Internal size (mm)	Store capacity	Weight (N.T/G.T)
HYC-950L	950 Liters	D780xW1200xH1894	D600xW1100xH1400	400 Bags	150Kg/180Kg



XC Series

Models: XC-88L/268L/358L/588L, Upright Type
Effective Capacity:
88/268/358/588 Liters
Temperature: 4°C



Temperature range 4±1°C • Vertical cabinet (upright) • Polyurethane insulation is 80-95 mm • Microprocessor-based temperature control system, printer is standard • Forced air circulation system, defrost automatic • Large screen LCD display of temperature, display accuracy of 0.1°C • Digital display of upper and lower temperature • Perfect alarm system: Audible/visual alarm, High or low temperature alarm, No battery alarm, Door ajar alarm, Filter blocking alarm, system failure • Exterior made from high quality steel plate, Inside is stainless steel • 3-layer transparent door with lock, door with heater around to prevent ice leak • Inner is drawer (stainless steel). 3units drawer for XC-88L, 4units for XC-268L, 5units for XC-358L, 6units for XC-588L • Interior fluorescent light • China good quality compressor (1 unit) and fan motor • The high-effective air-cooled condenser and the fin-type evaporator • Refrigerant as R134a, CFC free • Power supply: 220V/50Hz 1 phase, can change as 220V/60Hz or 110V/50Hz or 110V/60Hz • Certificate: CE mark, ISO9001, ISO14001, ISO13485 • Optional : chart recorder.

Model	Capacity	External size (mm)	Packing size (mm)	Blood capacity
XC-88L	88 Liters	D500xW450xH1515	D402xW343xH663	450ml/bag x 40-50/bags
XC-268L	268 Liters	D700xW628xH1662	D512xW518xH973	450ml/bag x 100-140/bags
XC-358L	358 Liters	D700xW628xH1962	D512xW515xH1270	450ml/bag x 160-200/bags
XC-588L	588 Liters	D800xW760xH1977	D550xW560xH1260	450ml/bag x 250-280/bags